



# Certainty can be effortless with Accuro

## Accuro®


### See what other ultrasound systems cannot.

Accuro's industry-first technology enables you to see what other ultrasound systems cannot see.

- **BoneEnhance®**, which is designed to visualize bone landmarks;
- and **SpineNav3D™**, which guides you with automatic epidural location and depth.

### Simply superior:

- Easy-to-use, automated image interpretation
- Superior imaging performance
- Portable, handheld form factor with rotatable touchscreen display
- Pocket size and battery operated

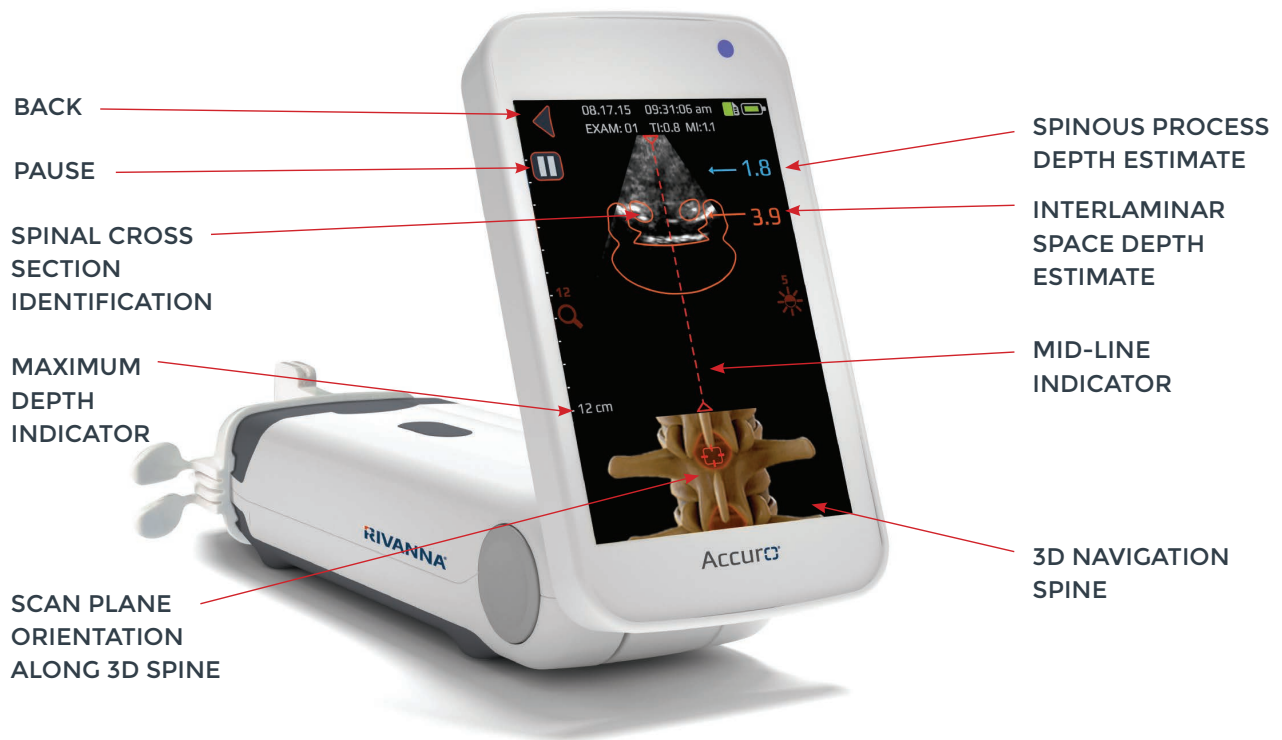



# Certainty can be effortless with Accuro

## Accuro®

**Intentionally Simple. No strings—or sonographer—attached.**

Accuro is a self-contained unit consisting of an ultrasound system, ultrasound probe, and rotatable touchscreen display.



### Specifications

#### #901-00002 Instrument specifications

ITEM	SPECIFICATION
Maximum dimensions	5.1 cm x 6.5 cm x 16.5 cm
Weight	0.7 lb (0.31 kg)
Display type	48.96 mm x 73.44 mm (480 pixels x 320 pixels, 166 dpi)

#### Battery specifications

ITEM	DESCRIPTION
Battery type	Lithium ion
Battery life	Up to 90 minutes of continuous imaging
Charge time	< 4 hours from an empty battery to full charge

#### Probe specifications

ITEM	DESCRIPTION
Diameter	12 mm
Frequency	5 MHz
Sector size	60
Max scan depth	12.5 cm

#### Battery charging components and accessories

USB wall charger
Micro-USB to USB cable
Acoustic coupling gel
Carrying case